




The Heavy Duty Anprolene option

The AN2000 features all metal construction, including a double walled design with a powder coated exterior and stainless steel interior. Wider door (14") is designed for larger procedure trays, coiled flexible scopes and bulky instruments.



CYCLE MONITORING

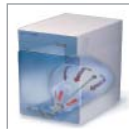

ACTIVE AERATION



FREE OPERATOR TRAINING



OPTIONAL ABATOR

Anprolene AN2000

The Professional Series still offers the most gentle **100% EtO sterilization** cycle on the market: **room temperature** sterilization (minimum 68° F), no steam injection, no deep vacuum, and **only 17.5 grams of gas**.

This is **the most robust and versatile low temperature sterilizer** available!

Anprolene sterilizers now features a number of important improvements. They include: more powerful ventilation and **purge pumps**, unique **active aeration** system, digital display **tracks cycle** progress, temperature sensor **monitors room conditions** during each cycle, microprocessor monitors cycle and automatically **restores system** in case of power failure.

AN2000 and Starter Kit



Suggested use:
Heavy duty use, wide trays or instruments, front office use.

Chamber dimensions:
20"L x 14"W x 14"H

Starter Kit - All you need for 40 cycles!

- 2 cases 20 cc ampoules and sterilization liner bags / **AN7916**
- 2 pouches Humidichip relative humidity stabilizer chips / **AN1071**
- 1 unit Humiditube holder / **AN1072**
- 2 boxes Dosimeter chemical indicators / **AN87**
- 2 boxes EtO exposure indicator strips - chemical indicators / **AN85**
- 2 boxes Steritest biological indicators / **AN80**
- 1 kit AirScan test kit (STEL and TWA) / **AN93**
- 1 unit Installation Kit / **AN74.01**

STERILIZER MUST BE PURCHASED WITH A STARTER KIT



Features and Specifications

Unmatched Sterilization Flexibility: The delicate Anprolene cycle protects and extends the life of your most delicate instruments; electronics, drills, fiber optics, cameras, laptop computers, and plastic or rubber items... Literally anything but food, drugs or liquids can be sterilized with Anprolene.

Sterilization and Aeration in the same cabinet: The Anprolene System utilizes an advanced purge system that aerates products for two hours at the end of the cycle. Active aeration continues indefinitely until products are removed, and the cycle is ended.

Unique Load Size Flexibility: The Anprolene system uses unique plastic sterilization liner bags and individual gas ampoules. We offer two different bag/ampoule combinations to match the size of your sterilization load. You never use more gas than you need.

Free Operator Training for the Life of the Sterilizer: We recommend that all Anprolene operators take our Key Operator Training course. Call to receive our study guide. When you are ready contact Andersen Customer Service to schedule a 20 minute test via telephone.

Proven Reliability: Anprolene is used in over 5,000 facilities worldwide. Find out how an Anprolene sterilizer can improve your practice today.

Specifications:

Model:
AN2000.14


Exterior dimensions:
22"L x 17.5"W x 22.5"H

Inner dimensions:
20"L x 14"W x 14"H

Total Weight:
74 pounds (33.6kg)

Power Supply:
115 +/-10%, 50/60 Hz/150 Watts

Manufactured by:
Andersen Sterilizers, Inc.

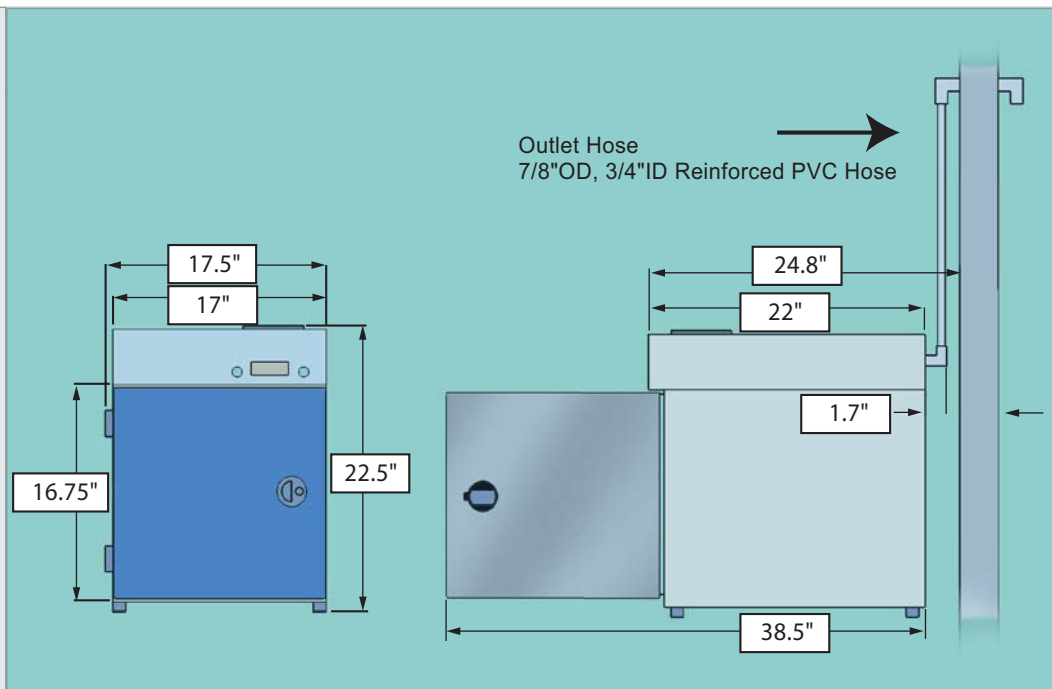

For more information about Anprolene® Sterilization Kits and prices call customer service at 800-523-1276, or visit www.anpro.com

